





Stand-alone and on-line PC operation

Direct parallel printer connection

Automatic test interpretation

Reusable or disposable flowmeter available

No calibration required



















## Spirobank Multifunction spirometer with alphanumeric display



MIR turbines are developed in full compliance with ATS standards and do not require calibration

## **Choose between two different flowmeters**

## Reusable turbine for long term operation

- High accuracy
- Long term stability
- Easy to clean or sterilize

## Disposable turbine for single patient use

- High accuracy
- Maximum hygiene guaranteed
- Easy replaceable without calibration
- Very low cost
- Available in boxes of 100 pieces





World's First!
Disposable turbine

International patent pending

## **Stand-alone spirometer**

- FVC, VC and MVV tests with PRE/POST comparison
- 26 parameters with automatic test interpretation and test quality control
- **Up to 100 test** memory capacity
- Internal temperature sensor for automatic BTPS conversion
- Several sets of predicted values
- Upgradeable internal software via internet at www.spirometry.com
- On-line PC connection

## **Easy telemedicine application**

• By using a special (optional) **internal software**, stored data can be transmitted via modem for **telemedicine applications** 

## **Direct printer connection**

- Stored test results can be printed by connecting the unit directly to a standard parallel printer.
- The **printout** gives a **full spirometry report** with **Flow/Volume curve**, results plus predicted values







# Spirobank Multifunction spirometer with graphic display


## **On-Line PC connection**

- WinspiroPRO software turns spirobank and spirobankG into an on-line clinical spirometer with the Flow/Volume curve shown on PC in real time
- PRE-POST bronchial challenge test
- **FEV1** dose-response curves
- Lung Age estimation
- Includes a series of user-selectable, and amusing, paediatric incentive animations



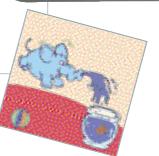

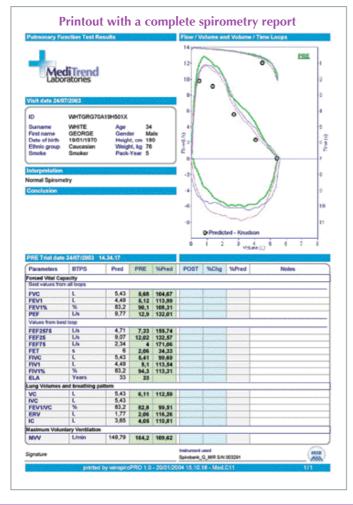

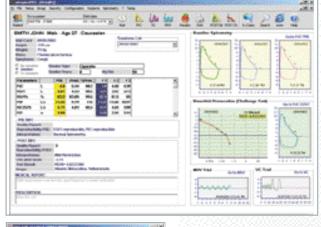



Note: spirobank and spirobankG have the same set of options, the only variation is the display

**Forced vital capacity:** FVC, FEV1, FEV1/FVC%, FEV6, FEV6%, PEF, FEF25, FEF50, FEF75, FEF25-75, FET, Vext, \*FVC, \*FEV1, \*PEF,

FIVC, FIV1, FIV1/FIVC%, PIF

**Slow vital capacity:** VC, IVC, IC, ERV, FEV1/VC% **Breathing pattern:** VT, VE, Rf, ti, te, ti/ttot, VT/ti

Max voluntary ventilation: MVV

(\* = Best value)

## **TECHNICAL SPECIFICATION**

**Temperature sensor:** semiconductor (0-45 °C)

Flow sensor: infrared interruption

**Max volume:** 10 L **Flow range:** ± 16 L/s

**Volume accuracy:** ± 3% or 50 mL, whichever is greater **Flow accuracy:** ± 5% or 200 mL/s, whichever is greater

**Dynamic resistance at 12L/s:** <0.5 cmH<sub>2</sub>O/L/s **Display (SpirobankG):** 120 x 32 pixel, graphic LCD

**Display (Spirobank):** 16 characters x 2 lines, alphanumeric LCD

Keyboard: 5 keys

Communication port: RS-232, bidirectional

**Power supply:** 9V DC (PP3 battery) **Dimensions:** 162 x 49 x 34 mm **Weight:** 180 grams (with battery)

## Standard equipment

- Reusable turbine for long term operation
- New winspiroPRO PC software
- RS-232 Interface cable
- 2 paper mouthpieces ø 30mm
- 1 plastic noseclip
- Carrying bag

## **Options**

- **Disposable turbine** for single patient use available in boxes of 100 pieces
- USB to RS232 converter cable
- RS232 to parallel printer converter
- External modem for telemedicine applications





### MIR

Via del Maggiolino, 125 00155 Roma - Italy tel. +39 06.22754777 fax +39 06.22754785 www.spirometry.com mir@spirometry.com